


REACTION APPARATUS

QR-4060 RESIN REACTION APPARATUS

Designed for pilot plant use. The apparatus is listed complete under each capacity. All covers are 3/4" thick plates with a 45/25 center neck, 3 29/26 side necks and a 19/22 side neck. The condenser is fully jacketed and has a 29/26 inner joint. The separatory funnel is graduated and has teflon stopcock and a 29/26 inner joint. The stirring rod has attached blades and is precision ground for a stirrer bearing with a 45/25 inner joint. The inlet tube is curved at the bottom and has a 29/26 inner joint. The clamp and stand assembly have adjustable feet.

	8500ML	12000ML	18000ML	33000ML	40000ML	50000ML
	7"x	8 1/4"x	10 1/4"x	12 1/2"x	12 1/2"x	12 1/2"x
<u>Description</u>	<u>17 1/2"</u>	<u>17 1/2"</u>	<u>17 1/2"</u>	<u>17 1/2"</u>	<u>23 1/2"</u>	<u>29 1/2"</u>
Complete	-110	-221	-332	-443	-554	-665
Flask Only	-111	-222	-333	-444	-555	-666
Cover Only	-112	-223	-334	-445	-556	-667
Funnel Only	-113	-224	-335	-446	-557	-668
Condenser Only	-114	-225	-336	-447	-558	-669
Inlet Adapter Only	-115	-226	-337	-448	-559	-670
Stirring Rod Only	-116	-227	-338	-449	-560	-671
Stirrer Bearing Only	-117	-228	-339	-450	-561	-672
Thermometer Only	-118	-229	-340	-451	-562	-673
Thermometer Adapter						
with Sleeve Only	-119	-230	-341	-452	-563	-674
Clamp & Stand						
Assembly Only	-120	-231	-342	-453	-564	-675

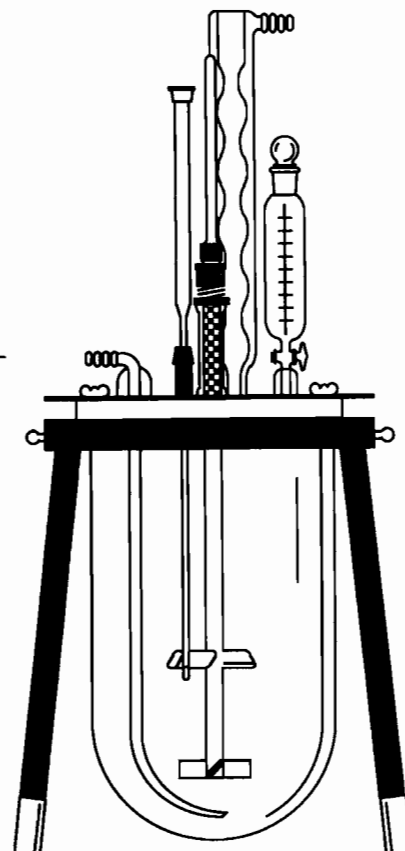

SEPTUMS

QS-3000 SEPTUM HOLDER

Tooled to hold a QS-3020. 7mm OD x 2.5mm ID tube approximately 6" long. Sold in packages of 12 each.

Code
-110



QS-3020 SEPTUM STOPPER

Red rubber sleeve style stopper, puncture seals automatically after needle is withdrawn. Sold in packages of 12 each.

<u>To Fit</u>	<u>Code</u>
QS-3000 Septum Holder or	
QS-3010 Septum Adapter	-110
14/20 or 14/35	-111
19/22 or 19/38	-112
24/25 or 24/40	-113

QS-3010 SEPTUM ADAPTER

For use when sampling air sensitive materials. Supplied with septum.

<u>Joint Size</u>	<u>Code</u>
14/20	-110
19/22	-111
24/40	-112

